



COOK MEDICAL EUROPE LTD.  
O'HALLORAN ROAD  
NATIONAL TECHNOLOGY PARK  
LIMERICK, V94 N8X2, IRELAND  
TEL: +353 61 334440 FAX: +353 61 334441  
WWW.COOKMEDICAL.EU

FSN & FSCA Ref: 2023FA0008

Date: 20JUL2023

## **Urgent Field Safety Notice – Medical Device Recall** **Lead Clippers**

For Attention of: Chief Executive / Risk Management / Purchasing

**Contact details of local representative (name, e-mail, telephone, address etc.)**

Cook Medical Europe Ltd.  
O'Halloran Road  
National Technology Park  
Limerick, Ireland  
E-mail: [European.FieldAction@CookMedical.com](mailto:European.FieldAction@CookMedical.com)  
Phone: Please refer to the attached Country Contacts List

For any further information or support concerning the information within this FSN,  
please contact your local Cook Medical Sales Representative or Cook Medical Europe  
Ltd.



COOK MEDICAL EUROPE LTD.  
O'HALLORAN ROAD  
NATIONAL TECHNOLOGY PARK  
LIMERICK, V94 N8X2, IRELAND  
TEL: +353 61 334440 FAX: +353 61 334441  
WWW.COOKMEDICAL.EU

FSN & FSCA Ref: 2023FA0008

## Urgent Field Safety Notice – Medical Device Recall Lead Clippers

### Risk addressed by FSN

<b>1. Information on Affected Devices</b>	
1.	<b>1. Device Type(s)</b> Lead Clippers are an auxiliary tool indicated for use in patients requiring percutaneous retrieval of cardiac leads.
1.	<b>2. Commercial name(s)</b> Lead Clippers
1.	<b>3. Primary clinical purpose of device(s)</b> Lead Clippers are used to separate the connector from pacemaker or defibrillator lead wire.
1.	<b>4. Device Model/Catalogue/Part Number(s)</b> Reference Part Number (RPN): LR-CLP001 Order Number (GPN): G20003
1.	<b>5. Affected serial or lot number range</b> N172912, N178633, N178677, N182691, N184187, N184320, N192312, N194720, N187302, N188171, N194884, N172737, N173178, N174134, N174864, N174906, N175201, N175277, N176323, N178518, N179397, N179728, N180115, N181397, N181437, N181629, N182095, N183533, N183826, N183939, N184517, N184881, N185570, N187777, N188290, N188322, N188490, N188741, N189965, N190539, N190720, N191357, N191583, N191973, N192341, N192497, N193021, N193571, N193731, N194362, N194666, N195113, N194084, N174170, N178093, N179534, N187519, N188615, N189751, N192120, N173653, N177898, N177797, N190400, N172875, N173959, N174841, N176830, N180459, N183042, N178609, N179301, N185880, N189099, N195038, N173977, N175964, N177515, N180077, N180953, N184990, N185239, N186318, N187000, N187727, N191791, N191870, N192310, N173691, N174294, N185335, N185615, N186683, N186963, N187107, N187359, N188466, N189666, N189805, N190813, N193490, N195447, N177102, N178814, N183261, N183457, N176716, N187246, N179021, N177239, N186888, N180504, N184518, N193119, N186749, and N189987



COOK MEDICAL EUROPE LTD.  
O'HALLORAN ROAD  
NATIONAL TECHNOLOGY PARK  
LIMERICK, V94 N8X2, IRELAND  
TEL: +353 61 334440 FAX: +353 61 334441  
WWW.COOKMEDICAL.EU

FSN & FSCA Ref: 2023FA0008

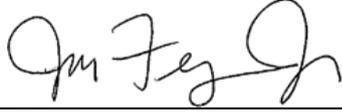
<b>2. Reason for Field Safety Corrective Action (FSCA)</b>	
<b>1. Description of the product problem</b>	
2.	Cook Medical Vandergrift identified that Lead Clippers may experience a complete breach of the chevron seal of the packaging. Therefore, the sterility of affected devices may be compromised.
<b>2. Hazard giving rise to the FSCA</b>	
2.	The affected devices may be non-sterile or contaminated with microorganisms. Potential adverse events that may occur if an affected product is used include infection, potentially being life-threatening and/or requiring medical/surgical intervention.  To date, Cook Medical Vandergrift has not received any customer complaints related to the adverse patient effects listed above for the affected lots. However, please be advised that compromised device sterility may go undetected by the user.

<b>3. Type of Action to Mitigate the Risk</b>	
<b>1. Action To Be Taken by the User</b>	
3.	<input checked="" type="checkbox"/> Identify Device <input checked="" type="checkbox"/> Quarantine Device <input checked="" type="checkbox"/> Return Device <input checked="" type="checkbox"/> Other Please complete the enclosed Customer Reply Form. Where product is indicated as being returned, our Customer Services department will contact you to organize the return and issue you with the relevant Returns Authorization number. Please include contact details on the Customer Reply form.  Returned Product should be addressed to: Cook Medical EUDC Robert-Koch-Straße, 2 52499 Baesweiler GERMANY  Credit will be provided for the returned affected products where applicable.
3.	<b>2. Is Customer Reply Required?</b> Yes. Form is attached specifying deadline for return.
3.	<b>3. Action Being Taken by the Manufacturer</b> <input checked="" type="checkbox"/> Product Removal
3.	<b>4. Is follow-up of patients or review of patients' previous results recommended?</b> Physicians should practice standard of care patient monitoring following the procedure for early identification of any complications to mitigate their severity. Cook Medical Vandergrift is not recommending additional patient monitoring as infection would likely present physical signs and symptoms abnormal to post-procedural patient recovery and promptly trigger medical intervention.



COOK MEDICAL EUROPE LTD.  
 O'HALLORAN ROAD  
 NATIONAL TECHNOLOGY PARK  
 LIMERICK, V94 N8X2, IRELAND  
 TEL: +353 61 334440 FAX: +353 61 334441  
 WWW.COOKMEDICAL.EU

FSN & FSCA Ref: 2023FA0008

4. General Information		
4.	1. FSN Type	New
4.	2. Further advice or information already expected in follow-up FSN?	No
4.	3. Manufacturer information For contact details of local representative refer to page 1 of this FSN.	
	a. Company Name	Cook Medical Vandergrift
	b. Address	1186 Montgomery Lane Vandergrift, PA 15690, United States
4.	4. The Competent (Regulatory) Authority of your country has been informed about this communication to customers.	
4.	5. Name/Signature	
		Jim Ferguson Quality Assurance Manager Cook Medical Vandergrift

Transmission of this Field Safety Notice
<p>This notice needs to be passed on all those who need to be aware within your organisation or to any organisation where the potentially affected devices have been transferred.</p> <p>Please transfer this notice to other organisations on which this action has an impact.</p> <p>Please maintain awareness on this notice and resulting action for an appropriate period to ensure effectiveness of the corrective action.</p> <p>Please report all device-related incidents to the manufacturer, distributor or local representative, and the national Competent Authority if appropriate, as this provides important feedback.</p>