



## McGRATH® MAC2 VIDEO LARYNGOSCOPE BLADE CLIP – LIMITED BATCH RECALL

4th November 2014

### BACKGROUND & SUMMARY

Aircraft Medical Limited have issued a voluntary recall on select batches of McGRATH® MAC2 Video Laryngoscope Blades due to an error in the manufacturing process that makes it possible for the blade retaining clip to be difficult to engage and thereby increase the possibility of a blade becoming disconnected during intubation, specifically when the device is being retracted from the patients mouth.

All affected blades should be quarantined immediately to avoid accidental use and shall be replaced free of charge by the Company with blades that are in compliance.

Only blades that display the recorded batch numbers are affected. Affected products should be returned to your distributor where they shall be decommissioned.

### AFFECTED PRODUCTS

Product Nomenclature :	McGRATH® MAC2 Video Laryngoscope Blade
Product Part Number :	350-017-000
Affected Products :	LOT 14030501
	LOT 14061616
	LOT 14031814
	LOT 13042401
	LOT 13122301
	LOT 13121202

### NATURE OF ISSUE

When a blade from an affected LOT is fitted to certain MAC Video Laryngoscopes it may be difficult to get the blade clip to engage onto the tooth on the handle. Where the problem does exist it is very apparent to the user.

### REQUIRED ACTION

All notified distributors shall check their own stock and their customers' stock for any blades that relate to the above LOT number, and arrange for any remaining stock to be removed from use. Note that it is possible that the blades have already been used without any problems.

#### Head Office

Aircraft Medical Limited  
10 Saint Andrew Square  
Edinburgh EH2 2AF  
United Kingdom  
Tel: +44 (0) 131 718 6042  
Fax: +44 (0) 131 718 6100  
www.aircraftmedical.com

Vat reg 801 6280 62  
Company SC 220767


#### Barcelona Office

Aircraft Medical (Barcelona), S.L.  
Llacuna 162 M 207  
08018 Barcelona  
Spain  
Tel: +34 93 309 06 68  
Fax: +34 93 401 97 62


CIF B-64202971

**IDENTIFICATION OF LOT NUMBER**

On McGRATH® MAC2 Peel Pouch

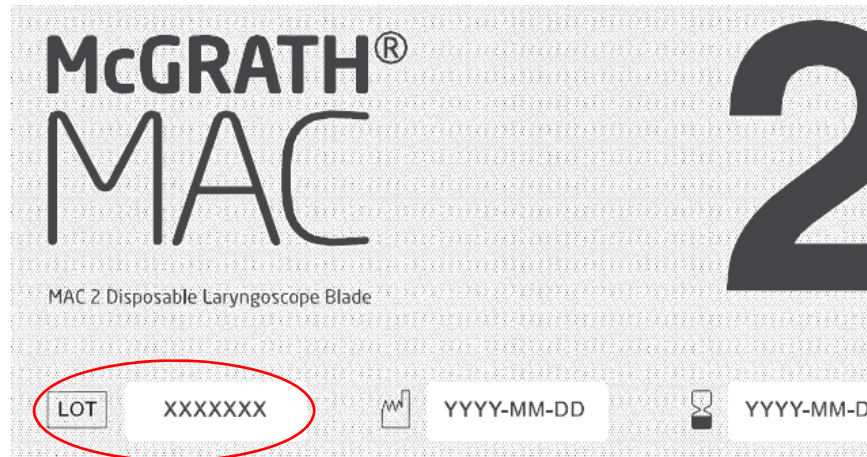

On McGRATH® MAC2 Box

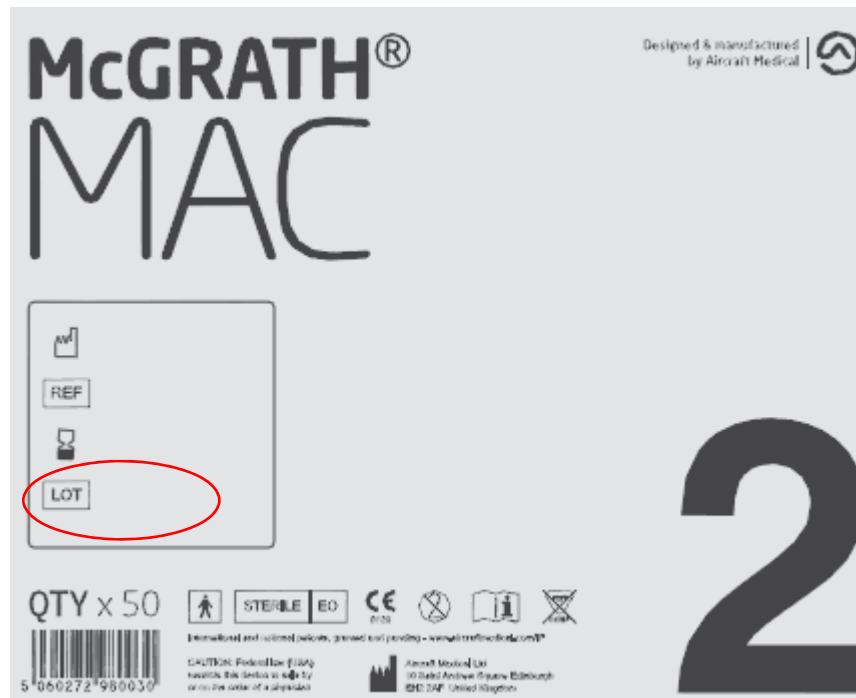

Thank you for your help and understanding on this matter.

**Jeremy Wheeler**  
Regulatory Affairs Director